

# Vcheck SDMA

The best choice for SDMA testing




# 1,

## Stronger correlation with the Gold standard method

The correlation of Vcheck SDMA ( $R^2=0.9624$ ) with HPLC is stronger than that of Analyzer 'C' of company 'I' ( $R^2=0.8871$ ).


\* HPLC (High-Performance Liquid Chromatography): a Gold standard of SDMA

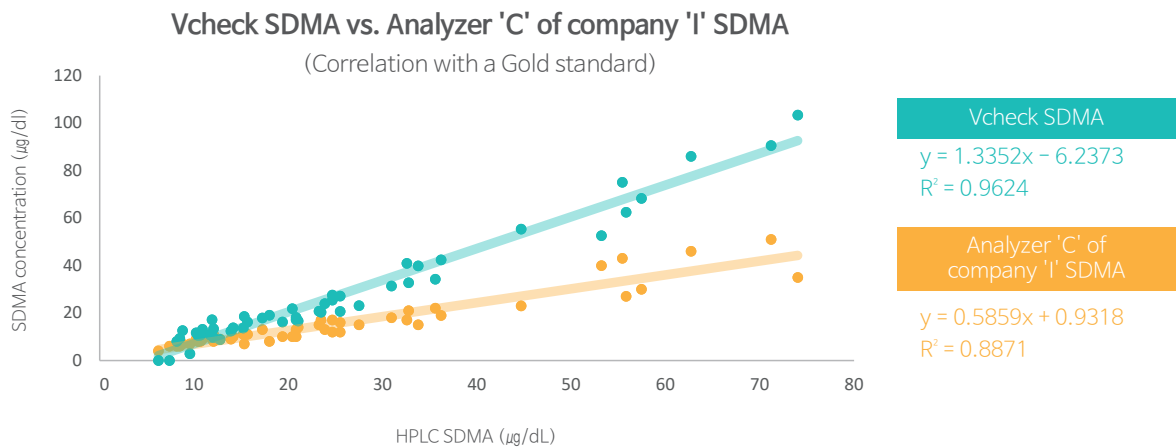

# 2,

## High correlation with company 'I' laboratory


With Vcheck SDMA, you can check the SDMA result values on the spot within 11 minutes. Vcheck SDMA has a high correlation (96.9%) with company 'I' laboratory.

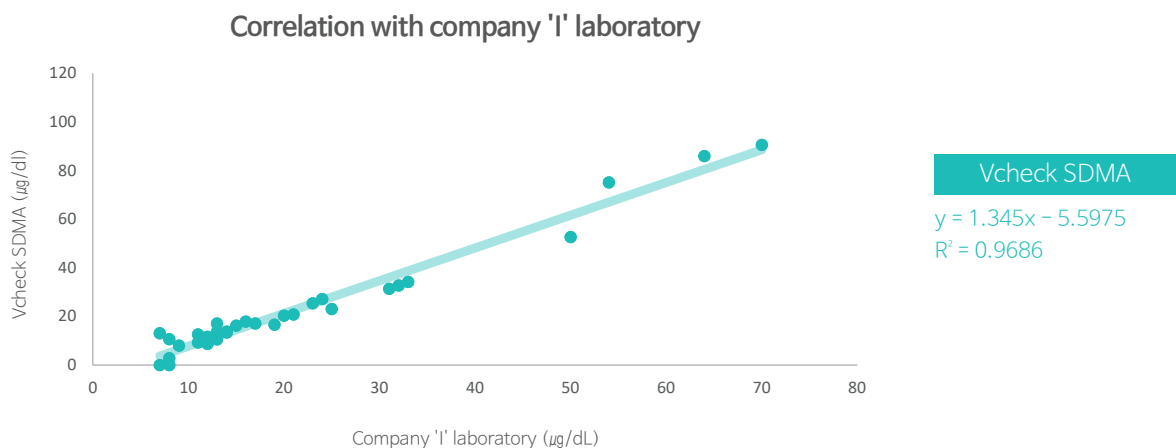

# 3,

## Excellent Spike-Recovery results


Spike-Recovery rate of Vcheck SDMA ( $R^2=0.9686$ ) is higher than that of analyzer 'C' of company 'I' ( $R^2=0.8793$ )

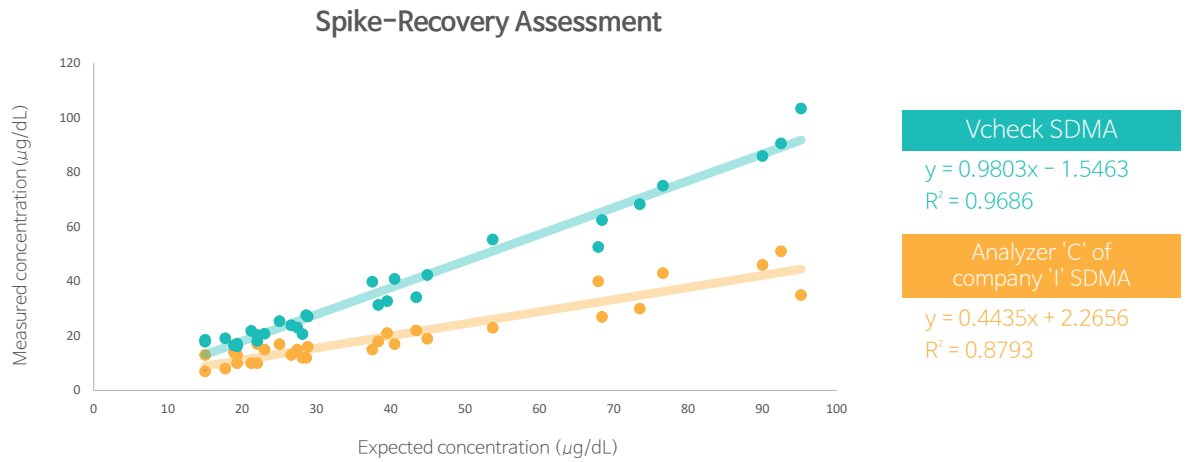

# 4,

## High specificity

Vcheck SDMA has been verified to have few cross-reactions with 3 analogues (NMMA, ADMA, L-arginine), which have similar structures to SDMA (less than 1%).

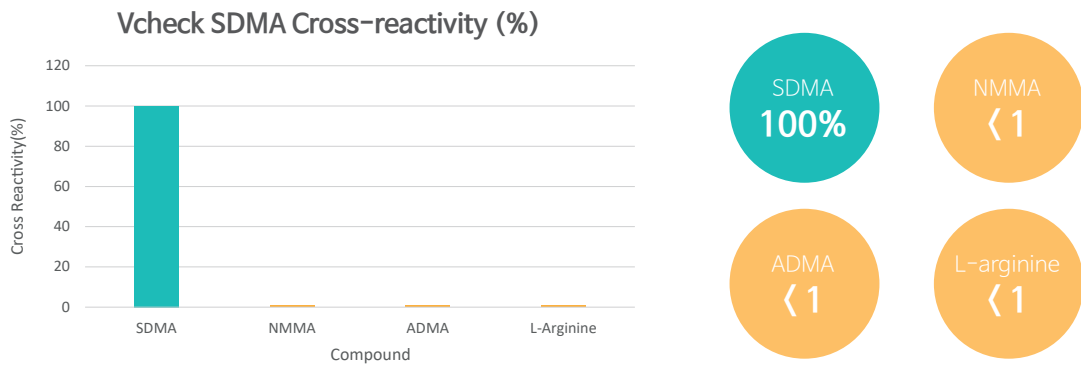

No need to wait for laboratory test results anymore!

**Cost effective and Accurate results,  
Vcheck SDMA makes it possible.**

