



Duration of antibody response following vaccination against feline immunodeficiency virus

Journal of Feline Medicine and Surgery
1–10

© The Author(s) 2016

Reprints and permissions:

sagepub.co.uk/journalsPermissions.nav

DOI: 10.1177/1098612X16673292

jfms.com

This paper was handled and processed by the European Editorial Office (ISFM) for publication in *JFMS*



Mark E Westman¹, Richard Malik², Evelyn Hall¹,
Matthew Harris³, Margaret J Hosie³ and Jacqueline M Norris¹

Abstract

Objectives Recently, two point-of-care (PoC) feline immunodeficiency virus (FIV) antibody test kits (Witness and Anigen Rapid) were reported as being able to differentiate FIV-vaccinated from FIV-infected cats at a single time point, irrespective of the gap between testing and last vaccination (0–7 years). The aim of the current study was to investigate systematically anti-FIV antibody production over time in response to the recommended primary FIV vaccination series.

Methods First, residual plasma from the original study was tested using a laboratory-based ELISA to determine whether negative results with PoC testing were due to reduced as opposed to absent antibodies to gp40. Second, a prospective study was performed using immunologically naive client-owned kittens and cats given a primary FIV vaccination series using a commercially available inactivated whole cell/inactivated whole virus vaccine (Fel-O-Vax FIV, three subcutaneous injections at 4 week intervals) and tested systematically (up to 11 times) over 6 months, using four commercially available PoC FIV antibody kits (SNAP FIV/FelV Combo [detects antibodies to p15/p24], Witness FeLV/FIV [gp40], Anigen Rapid FIV/FelV [p24/gp40] and VetScan FeLV/FIV Rapid [p24]).

Results The laboratory-based ELISA showed cats from the original study vaccinated within the previous 0–15 months had detectable levels of antibodies to gp40, despite testing negative with two kits that use gp40 as a capture antigen (Witness and Anigen Rapid kits). The prospective study showed that antibody testing with SNAP Combo and VetScan Rapid was positive in all cats 2 weeks after the second primary FIV vaccination, and remained positive for the duration of the study (12/12 and 10/12 cats positive, respectively). Antibody testing with Witness and Anigen Rapid was also positive in a high proportion of cats 2 weeks after the second primary FIV vaccination (8/12 and 7/12, respectively), but antibody levels declined below the level of detection in most cats (10/12) by 1 month after the third (final) primary FIV vaccination. All cats tested negative using Witness and Anigen Rapid 6 months after the third primary FIV vaccination.

Conclusions and relevance This study has shown that a primary course of FIV vaccination does not interfere with FIV antibody testing in cats using Witness and Anigen Rapid, provided primary vaccination has not occurred within the previous 6 months. Consequently, Witness and Anigen Rapid antibody test kits can be used reliably to determine FIV infection status at the time of annual booster FIV vaccination to help detect ‘vaccine breakthroughs’ and in cats that have not received a primary course of FIV vaccination within the preceding 6 months. The duration of antibody response following annual booster FIV vaccination and the resulting effect on antibody testing using PoC kits needs to be determined by further research. The mechanism(s) for the variation in FIV antibody test kit performance remains unclear.

Accepted: 9 September 2016

¹Sydney School of Veterinary Science, The University of Sydney, Sydney, NSW, Australia

²Centre for Continuing Veterinary Education, The University of Sydney, Sydney, NSW, Australia

³Centre for Virus Research, The University of Glasgow, Glasgow, UK

Corresponding author:

Mark E Westman BVSc, PhD, Sydney School of Veterinary Science, McMaster Building B14, The University of Sydney, Sydney, NSW 2006, Australia

Email: mark.westman@sydney.edu.au

Introduction

Feline immunodeficiency virus (FIV) can infect domestic cats and may cause, after a long asymptomatic phase, variable clinical disease due to its immunosuppressive and oncogenic properties.^{1,2} A FIV vaccine (Fel-O-Vax FIV; Boehringer Ingelheim), released for use in domestic cats (2002 in the USA; 2004 in Australia), was used as a 'proof of concept' for the development of sterilising immunity against lentiviruses, including human immunodeficiency virus 1.³⁻⁵ However, the vaccine was not registered in many jurisdictions (eg, Europe), owing, in part, to concerns related to the production of antibodies in FIV-vaccinated cats indistinguishable from those produced in response to natural FIV infection, such antibodies being the target for all point-of-care (PoC) test kits.⁶ Recently, we showed that some FIV antibody detection kits could differentiate FIV-vaccinated and FIV-infected cats under field conditions, reinforcing the complexity of antibody responses that occur following FIV vaccination.^{7,8}

The FIV genome is composed of approximately 9500 nucleotides, comprising three main open reading frames (*gag*, *pol* and *env*), encoding major capsid proteins (matrix, p15; capsid, p24; nucleocapsid, p7), viral enzymes (protease; reverse transcriptase; integrase) and envelope glycoproteins (transmembrane, gp40; surface, gp120), respectively.⁹ Areas of the genome capable of evoking host antibody response (B-cell epitopes) have been identified in the p15, p24, p7, gp40 and gp120 domains, with immunodominant epitopes located in the highly variable region (V3) of gp120.^{3,10} A cascade of antibody responses occurs following natural FIV infection, with antibodies to p24 and gp40 detectable, using Western blot, within 3 weeks of infection and antibodies to p15 detectable within 4 weeks of infection.^{8,11} Antibody production in FIV-infected cats persists for life, although antibody levels (particularly to p15 and p24) may wane in the terminal stage of infection.³ The complex nature of this antibody cascade has resulted in variable definitions of FIV positivity based on results from Western blotting, including: (i) presence of antibodies to gp120; (ii) antibodies to gp120 and at least one core protein (p7, p15 or p24); (iii) antibodies to at least two core proteins; or (iv) antibodies to three core proteins.^{11,12}

The only commercial FIV vaccine (Fel-O-Vax FIV) available consists of formalin-inactivated whole cells (IWC) and whole virus (IWV) suspended together in an adjuvant.⁴ It was presumed, on the basis of this composition, that the antibody response of FIV-vaccinated cats would be indistinguishable from that of FIV-infected cats.¹³ However, by using three PoC antibody tests detecting antibodies to different target FIV antigens, it was shown that FIV-vaccinated/FIV-uninfected cats consistently tested FIV positive with SNAP FIV/FelV Combo (p15 and p24) but FIV negative with Witness FelV/FIV (gp40) and Anigen Rapid FIV/FelV (p24 and gp40) for

0–7 years following FIV vaccination. Consequently, there was speculation that p15 retains immunogenicity during FIV vaccine production to a greater extent than p24 or gp40.⁷ This finding was questioned by another researcher, who used only Witness kits and a small cohort of kittens ($n = 19$) given a primary FIV vaccination series not in accordance with the manufacturer's recommendations (two injections administered instead of three), with a high FIV false-positive rate reported using Witness.¹⁴ Another larger study ($n = 104$), however, confirmed the ability of the Witness and Anigen Rapid kits to differentiate FIV-vaccinated and FIV-infected cats, but poor results from a fourth test kit that only detects antibodies to p24 (VetScan FeLV/FIV Rapid) challenged the notion that the ability of test kits to differentiate is solely linked to the choice of FIV antigen for antibody capture.¹⁵ Further work is therefore required to determine precisely and prospectively the antibody response following FIV vaccination in relation to PoC test kit methodology.

The aims of the current study were: (i) to determine if FIV-vaccinated cats produce antibodies to gp40 at concentrations below the detection threshold of Witness FelV/FIV and Anigen Rapid FIV/FelV kits, using a laboratory well-based ELISA; and (ii) to investigate semiquantitatively the duration of antibody response to p15, p24 and gp40 in cats following a primary course of FIV vaccination using four PoC FIV antibody test kits (SNAP FIV/FelV Combo, Witness FelV/FIV, Anigen Rapid FIV/FelV and VetScan FeLV/FIV Rapid).

Materials and methods

Sample population (study 1)

Residual blood from Westman et al was used for the first arm of this study.⁷ A total of 118 FIV-vaccinated cats had been recruited, comprising four FIV-infected and 114 FIV-uninfected cats. The median age of these cats was 7 years (range 2–18 years, interquartile range [IQR] 5–10 years) and the procedure for final assignment of FIV status, which included a combination of antibody testing, PCR testing and, occasionally, virus isolation, was described previously.⁷ All 118 cats had received a primary course of FIV vaccination consisting of three vaccines 2–4 weeks apart, in accordance with the manufacturer's recommendations, and a minimum of two annual booster vaccines, with no more than a 15 month gap between annual vaccinations. In total, 110/118 had received three or more annual boosters. Most cats (105/118) had been vaccinated within 1 year of sampling, and all cats had been vaccinated within 15 months of sampling (range 2–443 days, median 215 days, IQR 126–308 days) (Figure 1). Seven cats from the original study overdue for their annual FIV vaccination by 3–7 years were not tested.⁷ Of the 114 FIV-uninfected cats, 114 had tested FIV positive with SNAP Combo, six with Witness and none with Anigen Rapid. Each cat was only

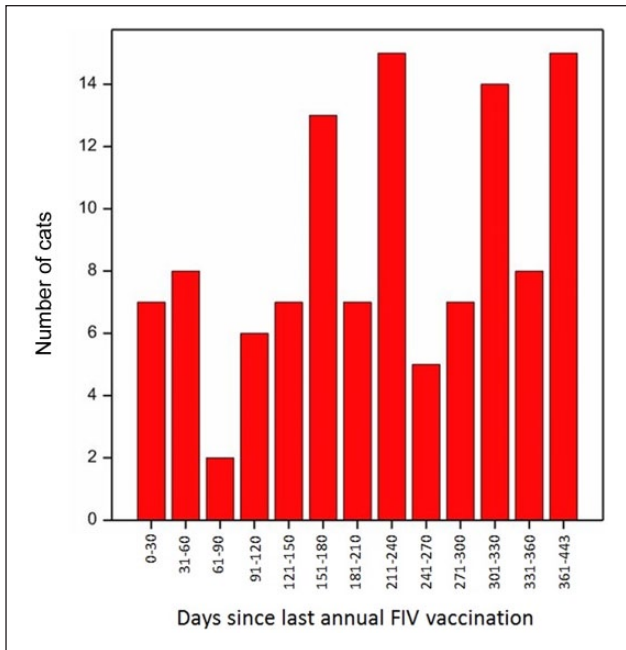


Figure 1 Categorisation of feline immunodeficiency virus (FIV) uninfected cats from study 1 based on time (days) elapsed since last FIV vaccination ($n = 114$). Of the 114 FIV-uninfected cats, 114 had tested FIV positive with SNAP Combo, six with Witness and none with Anigen Rapid. The six FIV false-positive results obtained with Witness occurred at the following intervals after FIV vaccination: 0–30 days ($n = 1$), 121–150 days ($n = 1$), 181–210 days ($n = 1$), 241–270 days ($n = 1$) and 331–360 days ($n = 2$)

available for sampling at a single time point, although occasionally a discordant cat had subsequent follow-up testing. A total of 23 FIV-unvaccinated/FIV-infected cats from the original study, determined by antibody and PCR testing and comprising a median age of 6 years (range 3–16 years, IQR 5–10 years), were also tested.⁷ Plasma stored at -80°C was transported on ice to Veterinary Diagnostic Services, The University of Glasgow for a laboratory-based gp40 ELISA. Approval was granted by The University of Sydney Animal Ethics Committee (approval number N00/1-2013/3/5920).

Sample population (study 2)

Four FIV-unvaccinated/FIV-uninfected kittens (<6 months of age) and 12 FIV-unvaccinated/FIV-uninfected cats (>6 months) were recruited from two veterinary clinics and two animal shelters in Sydney, Australia. The median age of all recruited cats was 2 years (range 0.3–8 years, IQR 1–4 years), significantly younger than cats in study 1 ($P < 0.001$; Mann–Whitney U-test). Recruited cats were given a primary course of three FIV vaccines subcutaneously, 4 weeks apart (weeks 0, 4 and 8), in accordance with the manufacturer’s recommendations, and antibody tested regularly (up to 11 times) using four

PoC FIV antibody test kits for 34 weeks (238 days; Table 1). Antibody testing at weeks 14, 16 and 20 was only pursued in cats that tested FIV positive with Witness or Anigen Rapid at the previous sampling, given the high likelihood of negative results with Witness/Anigen Rapid and positive results with SNAP Combo/VetScan Rapid in the other cats; one of these cats was lost to follow-up and unable to be tested at weeks 16 and 20. PCR testing was performed by a commercial laboratory (FIV RealPCR; IDEXX Laboratories) at the start of the study (week 0; prior to the first FIV vaccine being given), and at the end of the study (week 34), to ensure FIV infection had not occurred during the course of vaccinations and period of antibody testing.

Owners were offered free FIV testing and vaccination in return for enrolling their cat in the study. Cats were housed with their owners for the duration of the study; outdoor access was not regulated and was at the owners’ discretion. One cat tested FIV positive with an antibody test kit at week 0 (Anigen Rapid; this cat was inexplicably FIV negative on PCR testing at week 0 and week 34, and remained FIV positive with Anigen Rapid throughout the 34 weeks. Antibody testing at week 0 was negative with SNAP Combo, Witness and VetScan Rapid.) and was ultimately withdrawn at the conclusion of the study owing to uncertainty regarding its FIV status. Three other cats were withdrawn during the study for various reasons unrelated to blood sampling or FIV vaccination (one cat was hit by a car and died between week 0 and week 2; one was withdrawn at the owner’s request after week 4 owing to transport difficulties; and one cat was euthanased by the shelter after week 20 as the cat was re-surrendered following an incident of human-directed aggression at home). All cats tested feline leukaemia virus (FeLV) negative with the four kits. Approval was granted by The University of Sydney Animal Ethics Committee (approval number N00/1-2015/858).

Detection of antibodies to gp40 using a laboratory ELISA

A peptide ELISA, using a nine-amino acid sequence (CNQNQFFCK; cysteine–asparagine–glutamine–asparagine–glutamine–phenylalanine–phenylalanine–cysteine–lysine)³ from the highly conserved immunodominant TM2 domain of gp40, was used to detect antibodies.¹⁶ Plasma samples were first complement inactivated by incubation at 56°C for 30 mins. The wells of 96-well microtitre plates (Immulon 2 HB; Thermo Fisher Scientific) were coated with 250 ng/well of lyophilised gp40 epitope (AltaBioscience), diluted in sodium carbonate bicarbonate binding buffer (0.2 M anhydrous sodium carbonate, 0.2 M sodium carbonate and deionised water at a ratio of 1:11.5:4, respectively). The plates were incubated at 4°C overnight while being agitated at 30 rpm. The following day the wells were

Table 1 Outline of the prospective study design (study 2), showing time points for three primary feline immunodeficiency virus (FIV) vaccinations (weeks 0, 4 and 8), PCR testing (weeks 0 and 34) and antibody testing (up to 11 times between weeks 0 and 34). Antibody testing was not performed at weeks 14, 16 or 20 for cats that were FIV negative at week 12 with Witness and Anigen Rapid

	First FIV vaccine	2 weeks after first FIV vaccine	Second FIV vaccine	2 weeks after second FIV vaccine	Third (final) FIV vaccine	2 weeks after final FIV vaccine	4 weeks after final FIV vaccine	6 weeks after final FIV vaccine	8 weeks after final FIV vaccine	12 weeks after final FIV vaccine	26 weeks after final FIV vaccine
Time in weeks	0	2	4	6	8	10	12	14	16	20	34
Time in days	0	14	28	42	56	70	84	98	112	140	238
Diagnostic antibody testing	•	•	•	•	•	•	•	•	•	•	•
SNAP Combo	•	•	•	•	•	•	•	•	•	•	•
Witness	•	•	•	•	•	•	•	•	•	•	•
Anigen Rapid	•	•	•	•	•	•	•	•	•	•	•
VetScan Rapid	•	•	•	•	•	•	•	•	•	•	•
IDEXX RealPCR testing	•	•	•	•	•	•	•	•	•	•	•
FIV vaccination: three injections 4 weeks apart	•		•		•						

aspirated and washed five times with 200 µl phosphate buffered saline (PBS) supplemented with 0.1% Tween (PBST). Unabsorbed sites were blocked following incubation with 200 µl 2% low-fat milk powder in PBST (MP/PBST) for 1 h at room temperature. The wells were then aspirated and washed five times with 200 µl PBST, and 100 µl plasma added to the wells at a dilution of 1/200 (MP/PBST). The plates were sealed and incubated at room temperature for 1 h before being washed five times with 200 µl PBST, after which 100 µl biotinylated goat anticat secondary antibody (Vector Laboratories) was added to each well at a dilution of 1/1000 (MP/PBST). The plates were sealed and incubated at room temperature for 1 h. Wells were then aspirated and washed five times with 200 µl PBST, and 100 µl horseradish peroxidase conjugated to streptavidin added to each well at a dilution of 1/1000 (MP/PBST). The plates were sealed and incubated at room temperature for 20 mins, aspirated and washed five times with 200 µl PBST, and then 100 µl 3',3',5',5'-tetramethylbenzidine liquid (TMB Super Slow; Sigma Aldrich) added to each well. Plates were again sealed and incubated at room temperature for 30 mins before being read at 650 nm using a microplate reader (MultiSkan Ascent Plate Reader; MTX Lab Systems) and optical density (OD) values recorded.

Positive and negative controls were included on each test plate. The positive control plasma was collected from a cat, infected experimentally with the biological isolate of FIV_{GLS}, which tested FIV positive by Western blot and virus isolation. The negative control plasma was collected from an uninfected, specific-pathogen free cat that had been confirmed FIV negative by Western blot and virus isolation. ELISA results were not categorised as 'positive' or 'negative', but rather the range of antibody responses against gp40 were compared among the FIV-vaccinated cats tested.

Detection of antibodies using FIV PoC test kits

Blood was collected via jugular or cephalic venepuncture and stored in an EDTA tube at 4°C. FIV antibody testing was performed using four commercially available PoC kits within 24 h of sampling, in accordance with the manufacturers' recommendations. The kits tested were SNAP FIV/FeLV Combo (IDEXX Laboratories), Witness FeLV/FIV (Zoetis Animal Health), Anigen Rapid FIV/FeLV (BioNote) and VetScan FeLV/FIV Rapid (Abaxis). SNAP Combo is a lateral-flow ELISA kit, while the other three kits use immunochromatography to detect different FIV antibodies (Table 2). The fourth kit (VetScan Rapid) was added to the three kits tested in a previous study to include a methodology detecting antibodies to p24 alone.⁷ The results panel for each cat was photographed digitally at the time of testing. It is important to note that all four kits are marketed for the diagnosis of FIV infection, rather than the detection of antibodies produced in response to FIV vaccination as used in the current study.

Table 2 Feline immunodeficiency virus (FIV) target antigen for the antibodies detected using the four different point-of-care FIV antibody kits tested in study 2

FIV antibody detection kit	FIV target antigen		
	p15	p24	gp40
SNAP FIV/FelV Combo (Australia, New Zealand, North America)*	•	•	
Witness FeLV/FIV			•
Anigen Rapid FIV/FelV		•	•
VetScan FeLV/FIV Rapid		•	

*SNAP FIV/FelV Combo sold in Europe has an additional target antigen (gp40) included

Statistical analysis

Numerical analyses were performed at the conclusion of the study using statistical software (Genstat 16th Edition; VSN International). Significance was considered at $P < 0.05$. A Shapiro–Wilk test was used to assess data for normality; as data were not normally distributed (age of cats in study 1 and study 2, days post-FIV vaccination in study 1 and gp40 ELISA OD values) medians were reported and Mann–Whitney U-tests used for comparisons. ANOVA testing was used on \log_e transformed data to compare gp40 ELISA OD values grouped according to months since last annual FIV vaccination (0–3, 3–6, 6–9 and 9–15 months), number of annual booster FIV vaccinations administered (2–8) and age of cat at testing (grouped <5 years, 5–10 years, >10 years). Simple linear regression modelling was also performed with \log_e OD values as the outcome and days since last annual FIV vaccination, or age of cat at testing, as explanatory variables. Multivariate regression modelling was performed to consider the combined effect of days elapsed since last vaccination and number of annual booster vaccinations administered.

Results

FIV gp40 laboratory quantitative ELISA (study 1)

FIV-infected cats tested positive for antibodies recognising gp40, irrespective of FIV vaccination status ($P < 0.001$ compared with negative control). FIV-unvaccinated/FIV-infected cats showed a similar antibody response to the positive control ($P = 0.20$), as expected, while FIV-vaccinated/FIV-infected cats showed a weaker antibody response compared with the positive control ($P = 0.01$). No significant difference in the magnitude of antibody response was observed between FIV-unvaccinated/FIV-infected and FIV-vaccinated/FIV-infected cats ($P = 0.07$) (Figure 2a). FIV-vaccinated/FIV-uninfected cats that tested FIV true negative with Witness kits in the original study ($n = 108$) also tested antibody positive ($P < 0.001$ compared with negative control), but antibody levels for these cats were lower than in FIV-infected cats ($P < 0.001$). FIV-vaccinated/FIV-uninfected cats that tested FIV false

positive with Witness kits in the original study ($n = 6$) tested antibody positive ($P < 0.001$ compared with negative control), with higher antibody levels compared with the 108 Witness true negative cats ($P < 0.001$). This distinction, however, was not clear cut; for example, the upper range of the 108 Witness true negative cats encompassed the six Witness false-positive cats apart from one individual (Figure 2a). When time since last vaccination was analysed as a potential factor in the antibody response of the 114 FIV-vaccinated/FIV-uninfected cats, no significant effect was found ($P = 0.42$ [days], $P = 0.07$ [grouped by month]) (Figure 2b). When age of cat at testing was considered as a possible factor affecting the antibody response of the 114 FIV-vaccinated/FIV-uninfected cats, no significant effect was found ($P = 0.21$ [years], $P = 0.20$ [grouped by category <5, 5–10, >10 years]) (Figure 2c). There was no significant difference in antibody response when cats were grouped according to the number of annual FIV vaccinations administered ($P = 0.43$; Figure 2d). Similarly, when the time since the last vaccination and the number of annual vaccinations were considered together, no significant effect was observed ($P \geq 0.61$).

FIV PoC testing (study 2)

Sixteen cats commenced the study and were vaccinated against FIV, with 12/16 cats completing the study. Table 3 provides a summary of results for these 12 cats. The supplementary material provides a summary of results for all 16 cats (ie, including the four exclusions), results for kittens <6 months of age ($n = 4$) and results for cats >6 months ($n = 8$).

Considering the 12 cats, FIV antibodies were detected as early as 2 weeks after the first vaccination using SNAP Combo and Witness, and as early as 4 weeks using Anigen Rapid and VetScan Rapid. Two weeks after the second vaccination (week 6), all cats (12/12) tested FIV positive with SNAP Combo and VetScan Rapid, 8/12 (67%) tested FIV positive with Witness and 7/12 (58%) tested FIV positive with Anigen Rapid. At the completion of the study, 6 months after the third vaccination (week 34), all cats (12/12) were FIV positive with SNAP Combo, two cats had become FIV negative with VetScan Rapid and all cats were FIV negative with Witness and Anigen Rapid (Figure 3).

Three cats were tested between weeks 14 and 20 as a consequence of testing FIV positive with Witness and/or Anigen Rapid at week 12: (i) one cat tested FIV positive with Witness at weeks 12 and 14, then was lost to follow-up until week 34 when it tested FIV negative with Witness; (ii) one cat tested FIV positive with Witness and Anigen Rapid at weeks 12, 14 and 16 but FIV negative with both kits at week 20; and (iii) one cat tested FIV positive with Anigen Rapid at weeks 12 and 14 but FIV negative using this kit at week 20.

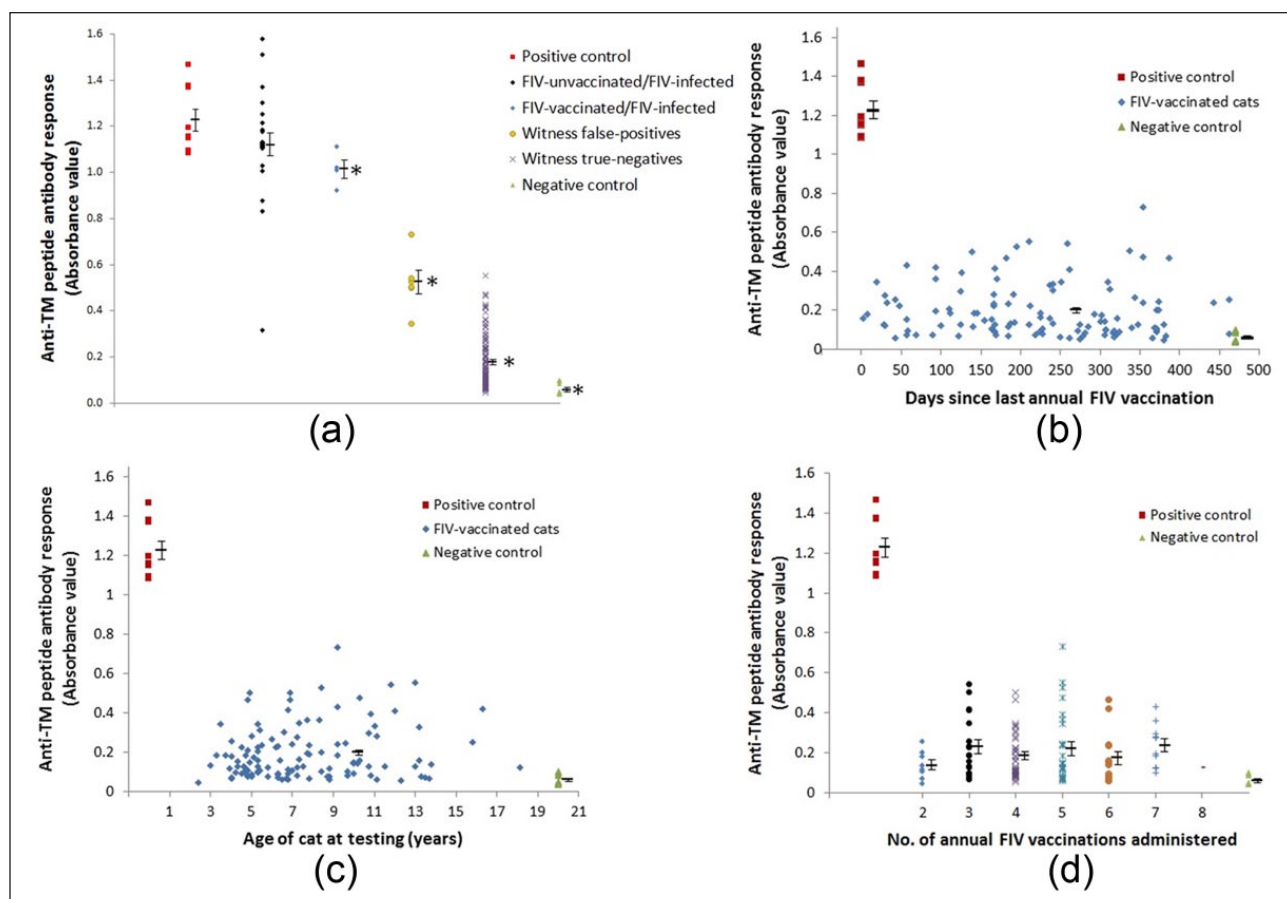


Figure 2 Results from ELISA testing for antibodies recognising feline immunodeficiency virus (FIV) gp40 peptide (study 1). Positive and negative controls are shown. The optical density (OD) is displayed on the y-axis. Mean and SEM bars are shown. (a) FIV-vaccinated/FIV-infected cats ($n = 4$), FIV-vaccinated/FIV-uninfected cats (FIV false positive with Witness, $n = 6$) and FIV-vaccinated/FIV-uninfected cats (FIV true negative with Witness, $n = 108$). *Significant difference ($P < 0.01$) between groups of cats. (b) FIV-vaccinated/FIV-uninfected cats ($n = 114$) according to the time (days) elapsed since the last annual FIV vaccination. No significant effect was found ($P = 0.42$). (c) FIV-vaccinated/FIV-uninfected cats ($n = 114$) according to the age of cat at the time of sampling (years). No significant effect was found ($P = 0.21$). (d) FIV-vaccinated/FIV-uninfected cats ($n = 114$) according to number of annual booster FIV vaccinations received. No significant effect was found ($P = 0.43$). TM = transmembrane

Discussion

The complexity of the antibody response following vaccination with a commercial IWC/IWV FIV vaccine was further described in this study. Laboratory-based ELISA quantification of antibodies to gp40 (study 1) demonstrated FIV-vaccinated cats had a detectable humoral response to gp40 for at least 15 months after FIV vaccination, despite a gp40 PoC test kit (Witness) testing negative in 95% of these samples (108/114).⁷ It was surprising not to find a quantitative decrease in gp40 antibody concentration over time since last FIV vaccination as determined by the ELISA OD value, especially when the Witness results from the second arm of the study (study 2) were considered. The explanation for this is unknown, and may relate to the older age of cats in study 1 compared with study 2, as well as reduced immunogenicity of the FIV vaccine with repeated booster vaccinations. Serial antibody testing using four different kits showed

that 6 months after a primary course of FIV vaccines was administered, Witness FeLV/FIV and Anigen Rapid FIV/FeLV tested FIV negative in 100% of cats, while SNAP FIV/FeLV Combo and VetScan FeLV/FIV Rapid tested FIV positive in 100% and 83% of cats, respectively.

At first glance, the Witness gp40 antibody results from study 2 and the original study appear contradictory.⁷ Study 2 found that a proportion of immunologically naive cats administered a primary course of FIV vaccination produced levels of antibodies to gp40 detectable by Witness for up to 6 months following vaccination, yet the original study found a very low level of FIV-positive results with Witness (6/114) in a cohort of cats vaccinated against FIV annually for at least 2 years, including an FIV-positive rate of only 1/16 in recently inoculated cats (vaccinated within the previous 12 weeks). Results obtained from the current study with Witness testing (study 2) were similar to results obtained by another

Table 3 Summary of feline immunodeficiency virus (FIV) PCR and FIV antibody test results from the prospective study (study 2) at various time points (n = 12). Four cats were withdrawn from the study (see text for details); results from these cats are not included here but are provided in the supplementary material

	First FIV vaccine	2 weeks after first FIV vaccine	Second FIV vaccine	2 weeks after second FIV vaccine	Third (final) FIV vaccine	2 weeks after final FIV vaccine	4 weeks after final FIV vaccine	6 weeks after final FIV vaccine	8 weeks after final FIV vaccine	12 weeks after final FIV vaccine	26 weeks after final FIV vaccine
Time in weeks	0	2	4	6	8	10	12	14	16	20	34
Time in days	0	14	28	42	56	70	84	98	112	140	238
Diagnostic antibody testing											
SNAP Combo	0/12	7/12	11/12	12/12	12/12	12/12	12/12	3/3	1/1	1/1	12/12
Witness	0/12	2/12	7/12	8/12	6/12	6/12	2/12	2/3	1/1	0/1	0/12
Anigen Rapid	0/12	0/12	1/12	7/12	5/12	4/12	2/12	1/3	1/1	0/1	0/12
VetScan Rapid	0/12	0/12	3/12	12/12	12/12	12/12	12/12	3/3	1/1	1/1	10/12
IIDEXX RealPCR testing	0/12										0/12

researcher who concluded that Witness testing alone could not be relied on to distinguish between natural FIV infection and FIV vaccination shortly after a primary FIV vaccination course.¹⁴ In that study, it was reported that 100% of FIV-vaccinated cats during their primary course tested FIV positive with Witness 4 weeks after the second FIV vaccination, 50% tested FIV positive 5 weeks after the second vaccination and 0% tested FIV positive 30 weeks after the second vaccination. There was a high dropout rate, however; of the 19 kittens that were enrolled, only four kittens were tested at 4 weeks, eight kittens were tested at 5 weeks and 11 kittens were tested at 30 weeks post-vaccination. Furthermore, only two FIV vaccinations were administered to kittens (instead of the recommended three), which made the results more difficult to interpret.¹⁴ A similar longitudinal study to the current design is required in adult cats prior to and following annual FIV vaccination to determine whether this period of detectable antibody response with PoC test kits such as Witness extends beyond primary FIV vaccination.

The explanation for this seeming discrepancy with Witness testing is possibly two-fold. First, some studies have reported a lower antibody response in people being re-vaccinated compared with those being vaccinated for the first time. Govaert et al found older people (>60 years) re-vaccinated with an inactivated influenza vaccine had a 'strikingly' lower humoral immune response than people who had not previously been vaccinated.¹⁷ A longitudinal study of elderly people (56–79 years) administered a primary pneumococcal polysaccharide vaccine, and given a booster vaccination 6 years later, found the antibody titres after re-vaccination were about half the titres after primary vaccination.¹⁸ The trend for a weaker antibody response following booster vaccination compared with primary vaccination, however, is not steadfast; for example, one study that investigated the vaccine-induced antibody response against hepatitis virus B in humans reported higher antibody levels following booster vaccination than at the end of the primary vaccination course 6 years earlier.¹⁹ Factors related to the nature of the pathogen and the antigen(s) and adjuvant presented during vaccination are likely to play a crucial role in determining the ongoing humoral response. Second, age may be a factor; cats recruited for the original Westman et al study (which became the cats in study 1)⁷ were substantially older than cats in study 2 (median age 7 vs 2 years; $P < 0.001$) owing to the large proportion of kittens (4/12). Likewise, Lappin tested only kittens.¹⁴ Duration of immunity (DOI) studies are sparse in the veterinary literature, and most are only concerned with protection from challenge rather than antibody quantitation for diagnostic purposes.²⁰ Flow cytometry studies have demonstrated an age-related remodelling of the immune system in cats, with a gradual decline in relative percentage of lymphocytes,²¹ and an absolute reduction in B cells in senior cats (10–14 years) compared

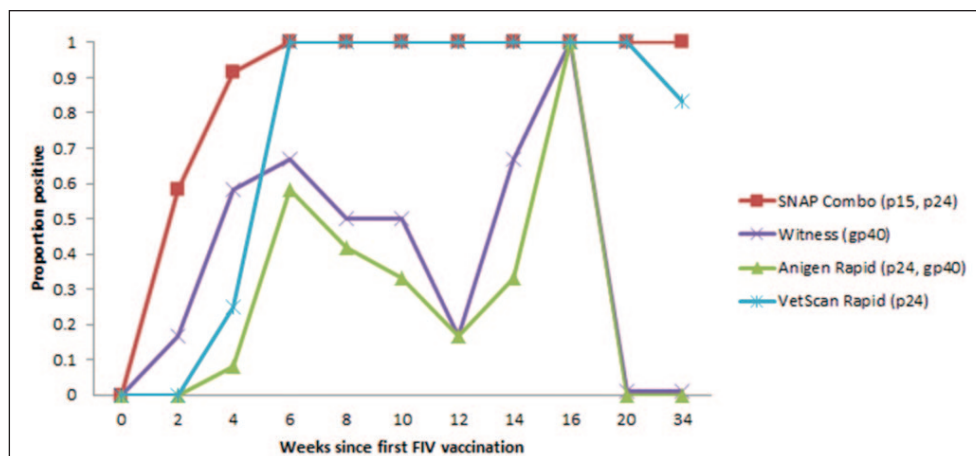


Figure 3 Summary of feline immunodeficiency virus (FIV) antibody test results from the prospective study (study 2) at various time points ($n = 12$). The FIV target capture antigen(s) for each point-of-care antibody test kit is included in brackets. A primary FIV vaccination course was administered at 0, 4 and 8 weeks

with young cats (2–5 years).²² DOI studies are more common in the human literature, where it is generally accepted that older people have a weaker humoral response following vaccination than younger people. For example, one study investigating antibody response in humans administered an inactivated H1N1 vaccine found pre- and post-vaccination titres were generally lower in the elderly (>70 years of age) than the young (<30 years).²³ We postulate that the accuracy of the Witness kit to assign correctly FIV infection status in FIV-vaccinated cats reported previously may be explained by a relatively low level of gp40 antibodies in older cats following booster FIV vaccination, compared with the younger cats in both study 2 and the study reported by Lappin,¹⁴ which were vaccinated against FIV for the first time. Contrary to this theory is the absence of a trend in study 1 for gp40 antibody concentration to diminish with increased age of cat at testing and/or number of annual FIV booster vaccinations administered (Figure 2c,d), nor was there a noticeable trend for kittens to test false positive with Witness more often than adult cats in study 2 (supplementary material). Inadequate sample sizes for both studies may have been responsible. To further investigate the role of re-vaccination and age, we plan to monitor cats in study 2 for several years and determine their gp40 antibody response following booster FIV vaccination, to see if their antibody response following annual vaccination is less than their antibody response following initial (primary) vaccination.

Sequential semiquantitative antibody testing following FIV vaccination with Witness and Anigen Rapid showed peak antibody production occurred during and shortly after a primary course of FIV vaccination (three injections at 4 week intervals). Two weeks after the second vaccination (week 6), 67% (8/12) and 58% (7/12) of cats tested seropositive for FIV antibodies using Witness and Anigen Rapid, respectively. By 4 weeks after the third vaccination (week 12), only 17% (2/12) of cats tested

seropositive with Witness or Anigen Rapid, and by 6 months after the third vaccination (week 34) p24 and gp40 antibody levels had decreased below the detection limit for both kits. Peak antibody production to p24 and gp40 6–12 weeks after the first primary FIV vaccination (in a series of three), as demonstrated by these results, supports results from previous studies, including results from p24 ELISA determinations in experimental kittens vaccinated three times 2–3 weeks apart (unpublished data, Boehringer Ingelheim). Using laboratory-based ELISA testing, Huang et al showed antibody against p24 and gp40 peaked 1–3 weeks after the third primary FIV vaccination (vaccines given 3 weeks apart) in kittens and decreased over the following 3 months, with antibody levels maintained for 12 months in most cats.^{24,25} The same result was found by another group using ELISA testing for antibodies to whole FIV antigen and recombinant p24 (r-gag) in cats aged 7–12 months.²⁶ Western blot analysis of four FIV-vaccinated cats confirmed antibody production to p15, p24 and gp40 3 weeks after the second primary FIV vaccination, which persisted for at least 12 months following the third primary FIV vaccination.¹³ The reason why ELISA gp40 testing in the current study (study 1) did not show a peak (and subsequent fall) in antibody production according to time elapsed since last vaccination (Figure 2b), similar to ELISA testing by Huang et al and Kusuhara et al,^{24–26} is uncertain. It may also relate to a reduced antibody response in older cats following annual booster FIV vaccination, rather than younger cats receiving a primary course of FIV vaccination.^{17,23}

Our results confirmed that care needs to be exercised in the period immediately following primary FIV vaccination using Witness and Anigen Rapid, with false-positive FIV results occurring using both. In the light of these findings, we suggest an amendment to our previous conclusion,⁷ and recommend that antibody testing

to detect FIV infection in FIV-vaccinated cats is reliable using Witness and Anigen Rapid, providing primary vaccination against FIV has not occurred within the preceding 6 months. As discussed above, more research needs to be performed to investigate whether this recommendation also applies to cats receiving annual booster FIV vaccinations. In a shelter situation, where large-scale FIV screening is being undertaken, it is highly unlikely cats will have received a primary FIV vaccination series in the preceding 6 months as FIV vaccination rates are generally low, and cats so vaccinated are generally well cared for and less likely to be surrendered to a shelter facility.²⁷ Transient antibody production for up to 6 months after primary FIV vaccination is also not relevant when testing cats for FIV 'vaccine breakthrough' immediately prior to the next annual FIV vaccination booster. If a positive FIV antibody result is obtained in a cat where recent primary FIV vaccination is possible, submitting seropositive specimens for confirmatory FIV PCR testing is recommended. A negative FIV test result with either Witness or Anigen Rapid remains robustly reliable and is recommended as the screening test of choice, except in cases of recent infection, when repeat testing 2 months later is recommended.²⁸

Results from sequential antibody testing in this study challenge our previous notion that p15 is more immunogenic than p24 and gp40 in the FIV vaccine.⁷ SNAP Combo (which detects antibodies to both p15 and p24) gave a seropositive result in 12/12 vaccinated cats from 6 weeks after the first FIV vaccination and all 12 cats remained FIV positive for the duration of the study (34 weeks). Additionally, VetScan Rapid (which detects antibodies to p24 but not p15) tested seropositive in 12/12 vaccinated cats from 6 weeks after the first FIV vaccination and 10/12 (83%) remained FIV positive at the end of the study. If the difference in performance between SNAP Combo and Witness/Anigen Rapid was solely attributable to p15 being more immunogenic in the FIV vaccine, then VetScan Rapid would have performed comparably well to Witness/Anigen Rapid in the current study. The differing performance of VetScan Rapid in study 2 infers that the difference in kit performance may rely more on factors related to testing methodology (eg, ELISA vs immunochromatography, the antibody threshold at which the test is set and how the capture antigen is prepared) than factors related to the FIV vaccine (ie, immunogenicity of different epitopes). For this reason, care must be taken when selecting an antibody kit to avoid false-positive results in FIV-vaccinated cats. Our findings cannot be extrapolated to other antibody kits without appropriate testing being performed. Future research should ideally quantitate the antibody response directed against each epitope (p15, p24 and gp40) over time, using a common methodology for each (eg, ELISA testing), to further understand the

breadth, magnitude and duration of the antibody cascade following FIV vaccination.

Conclusions

The complexity of antibody production following FIV vaccination was further described using both laboratory-based ELISA and an extended range of PoC test kits. Antibodies to p15, p24 and gp40 were detectable early (within 4 weeks of the first FIV vaccination) using various test kits. SNAP Combo and VetScan Rapid tested persistently FIV positive for 6 months in cats given a primary course of FIV vaccination, while Witness and Anigen Rapid tested FIV negative in all cats by 6 months following primary FIV vaccination. The limit of detection at which these antibody kits are calibrated appears to be the critical factor, as antibodies to gp40 (and likely p15 and p24) persist for at least 15 months after FIV vaccination and kits that are biased towards sensitivity (eg, SNAP Combo, VetScan Rapid) will detect these antibodies in addition to those produced by natural FIV infection. In jurisdictions where FIV vaccination is practised, testing for potential FIV breakthrough infection prior to annual FIV vaccination is prudent.²⁹ In this setting, Witness and Anigen Rapid are the FIV antibody test kits of choice as fewer false-positive results would be anticipated than with the SNAP Combo and VetScan Rapid test kits.

Acknowledgements We are grateful to the significant in-kind contribution of diagnostic kits by IDEXX Laboratories, Zoetis Animal Health, BioNote and Abaxis. Mark Westman was supported by an Australian Postgraduate Award (APA), the Neil and Allie Lesue Scholarship and the Herbert Johnson Travel Grant Scholarship. Richard Malik was supported by the Valentine Charlton Bequest. Alan Marcus and Phillip McDonagh (Head of Regulatory Affairs for Animal Health, Boehringer Ingelheim Australia) provided helpful comments to the manuscript. Finally, we thank all of the cats and owners who supported this study.

Conflicts of interest The authors declared no potential conflicts of interest with respect to the research, authorship, and/or publication of this article.

Funding This study was supported financially by the Australian Companion Animal Health Foundation (ACAHF), Feline Health Research Fund (FHRF) and Boehringer Ingelheim, Australia.

Supplementary Material The following files are available: Summary of FIV PCR and FIV antibody test results from the prospective study (study 2) at various time points (n = 16), including four cats withdrawn from the study at various time points (after weeks 0, 4, 20 and 34). Summary of FIV PCR and FIV antibody test results from the prospective study (study 2) at various time points for kittens <6 months of age (n = 4). Summary of FIV PCR and FIV antibody test results from the prospective study (study 2) at various time points for cats >6 months of age (n = 8).

References

- 1 Pedersen N, Ho E, Brown M, et al. **Isolation of a T-lymphotropic virus from domestic cats with an immunodeficiency-like syndrome.** *Science* 1987; 235: 790–793.
- 2 Shelton GH, Grant CK, Cotter SM, et al. **Feline immunodeficiency virus and feline leukemia virus infections and their relationships to lymphoid malignancies in cats: a retrospective study (1968–1988).** *J Acquir Immune Defic Syndr* 1990; 3: 623–630.
- 3 Lecollinet S and Richardson J. **Vaccination against the feline immunodeficiency virus: the road not taken.** *Comp Immun Microbiol Infect Dis* 2008; 31: 167–190.
- 4 Yamamoto JK, Pu RY and Sato E. **Feline immunodeficiency virus pathogenesis and development of a dual-subtype feline immunodeficiency-virus vaccine.** *AIDS* 2007; 21: 547–563.
- 5 Bienze D. **FIV in cats – a useful model of HIV in people?** *Vet Immunol Immunopathol* 2014; 159: 171–179.
- 6 Hosie MJ, Addie D, Belak S, et al. **Feline immunodeficiency: ABCD guidelines on prevention and management.** *J Feline Med Surg* 2009; 11: 575–584.
- 7 Westman ME, Malik R, Hall E, et al. **Determining the feline immunodeficiency virus (FIV) status of FIV-vaccinated cats using point-of-care antibody kits.** *Comp Immun Microbiol Infect Dis* 2015; 42: 43–52.
- 8 Yamamoto JK, Sparger E, Ho EW, et al. **Pathogenesis of experimentally induced feline immunodeficiency virus infection in cats.** *Am J Vet Res* 1988; 49: 1246–1258.
- 9 Olmsted RA, Hirsch VM, Purcell RH, et al. **Nucleotide sequence analysis of feline immunodeficiency virus – genome organisation and relationship to other lentiviruses.** *Proc Natl Acad Sci U S A* 1989; 86: 8088–8092.
- 10 Yamamoto JK, Sanou MP, Abbott JR, et al. **Feline immunodeficiency virus model for designing HIV/AIDS vaccines.** *Curr HIV Res* 2010; 8: 14–25.
- 11 O'Connor TP, Tanguay S, Steinman R, et al. **Development and evaluation of immunoassay for detection of antibodies to the feline T-lymphotropic lentivirus (feline immunodeficiency virus).** *J Clin Microbiol* 1989; 27: 474–479.
- 12 Hosie MJ and Jarrett O. **Serological responses of cats to feline immunodeficiency virus.** *AIDS* 1990; 4: 215–220.
- 13 Uhl EW, Heaton-Jones TG, Pu R, et al. **FIV vaccine development and its importance to veterinary and human medicine: FIV vaccine 2002 update and review.** *Vet Immunol Immunopathol* 2002; 90: 113–132.
- 14 Lappin MR. **Detection of feline immunodeficiency virus antibodies in serum of vaccinated cats using a commercially available kit (2015 ACVIM Forum Research Abstract Program).** *J Vet Intern Med* 2015; 29: 1122–1256.
- 15 Crawford C. **Does a DIVA test exist for differentiating FIV infection from FIV vaccination? (2016 ACVIM Form Research Abstract Program).** *J Vet Intern Med* 2016; 30: 1407–1519.
- 16 Avrameas A, Strosberg AD, Moraillon A, et al. **Serological diagnosis of feline immunodeficiency virus infection based on synthetic peptides from *env* glycoproteins.** *Res Virol* 1993; 144: 209–218.
- 17 Govaert TME, Sprenger MJW, Dinant GJ, et al. **Immune response to influenza vaccination of elderly people. A randomized double-blind placebo-controlled trial.** *Vaccine* 1994; 12: 1185–1189.
- 18 Mufson MA, Hughey DF, Turner CE, et al. **Revaccination with pneumococcal vaccine of elderly persons 6 years after primary vaccination.** *Vaccine* 1991; 9: 403–407.
- 19 Trivello R, Chiaramonte M, Ngatchu T, et al. **Persistence of anti-HBs antibodies in health care personnel vaccinated with plasma-derived hepatitis B vaccine and response to recombinant DNA HB booster vaccine.** *Vaccine* 1995; 13: 139–141.
- 20 Schultz RD. **Duration of immunity for canine and feline vaccines: a review.** *Vet Microbiol* 2006; 117: 75–79.
- 21 Heaton PR, Blount DG, Mann SJ, et al. **Assessing age-related changes in peripheral blood leukocyte phenotypes in domestic shorthaired cats using flow cytometry.** *J Nutr* 2002; 132: 1607S–1609S.
- 22 Campbell DJ, Rawlings JM, Koelsch S, et al. **Age-related differences in parameters of feline immune status.** *Vet Immunol Immunopathol* 2004; 100: 73–80.
- 23 de Bruijn IA, Remarque EJ, Jol-van der Zijde CM, et al. **Quality and quantity of the humoral immune response in healthy elderly and young subjects after annually repeated influenza vaccination.** *J Infect Dis* 1999; 179: 31–36.
- 24 Huang C, Conlee D, Loop J, et al. **Efficacy and safety of a feline immunodeficiency virus vaccine.** *Anim Health Res Rev* 2004; 5: 295–300.
- 25 Huang C, Conlee D, Gill M, et al. **Dual-subtype feline immunodeficiency virus vaccine provides 12 months of protective immunity against heterologous challenge.** *J Feline Med Surg* 2010; 12: 451–457.
- 26 Kusuhara H, Hohdatsu T, Okumura M, et al. **Dual-subtype vaccine (Fel-O-Vax FIV) protects cats against contact challenge with heterologous subtype BFIV infected cats.** *Vet Microbiol* 2005; 108: 155–165.
- 27 Salman MD, New JJG, Scarlett JM, et al. **Human and animal factors related to relinquishment of dogs and cats in 12 selected animal shelters in the United States.** *J Appl Anim Welf Sci* 1998; 1: 207–226.
- 28 Levy J, Crawford C, Hartmann K, et al. **2008 American Association of Feline Practitioners' feline retrovirus management guidelines.** *J Feline Med Surg* 2008; 10: 300–316.
- 29 Westman ME, Malik R, Hall E, et al. **The protective rate of the feline immunodeficiency virus vaccine: an Australian field study.** *Vaccine* 2016; 34: 4752–4758.